ThevoSchlummerstern Senso

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Dear Customer,

Thank you for choosing our product and for the trust you have placed in us. In order for you and the user to enjoy your new product safely, conveniently, and comfortably, we kindly ask that you read the Instructions for Use. Doing this will also help avoid damage, injury, and incorrect operation and handling of the product. Please ensure that a copy of the Instructions for Use remains with the product so that every user has access to it.

If you have any questions, please contact your authorised distributor or us directly.

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Phone: +49 (0) 4761 8860

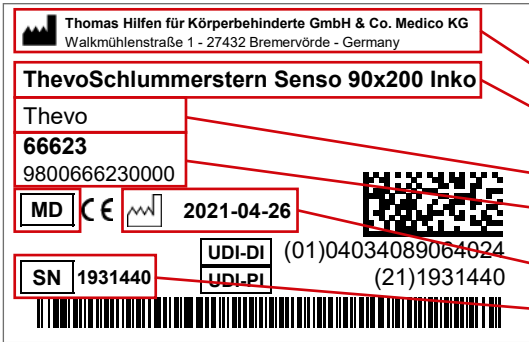
www.thomashilfen.de

We reserve the right to make technical changes. This can lead to differences between your product and how it is presented in brochures and other media.

Table of Contents

Figures	3
Preface / Manufacturer	12
Labels / Model Identification	13
General Remarks / Safety Instructions	14
Intended Purpose	14
– Intended Use / Application Areas	14
– Indication / Patient Target Group / Intended Users	14
– Risks Associated with Use and Contraindications	15
Product Description	15
– Scope of Delivery / Standard Equipment	15
– Options	15
– Putting into Service	15
– Teaching the User	15
Adjustment Options / Features	15
– Accessories	15
– Technical Data	17
– Material Specifications	17
Cleaning, Maintenance, and Repair	17
– Cleaning, Care Instructions, and Disinfection on a Daily Basis	17
– Maintenance and inspection by the distributor / Maintenance Schedule	18
– Spare Parts	18
– Storage / Disposal	18
EU Declaration of Conformity	18
Warranty	18
Period of Use / Service Life	19
Reuse	19
Distributor	19

Labels / Model Identification

Label	Meaning
 <p>Thomas Hilfen für Körperbehinderte GmbH & Co. Medico KG Walkmühlenstraße 1 - 27432 Bremervörde - Germany</p> <p>ThevoSchlummerstern Senso 90x200 Inko</p> <p>Thevo</p> <p>66623</p> <p>9800666230000</p> <p>MD CE 2021-04-26</p> <p>UDI-DI (01)04034089064024</p> <p>UDI-PI (21)1931440</p> <p>SN 1931440</p>	<p>Barcode Label</p> <p>Position on the Product (Fig. 1 + 2)</p> <ul style="list-style-type: none"> - Manufacturer - Product Name - Product Group - Item Number short / long - MD = Medical Device - Manufacture Date - Serial Number
 <p>Thomashilfen ThevoSchlummerst. S Kopfende hier / head end here</p> <p>Service</p> <p>www.thomashilfen.de</p> <p>Thomashilfen ThevoSchlummerst. S < 60 kg</p> <p>Service</p> <p>www.thomashilfen.de</p> <p>Thomashilfen für Körperbehinderte GmbH & Co. Medico KG Walkmühlenstraße 1 - 27432 Bremervörde - Germany</p> <p>ThevoSchlummerstern Senso 90x200 Inko</p> <p>Thevo</p> <p>66623</p> <p>9800666230000</p> <p>MD CE 2021-04-26</p> <p>UDI-DI (01)04034089064024</p> <p>UDI-PI (21)1931440</p> <p>SN 1931440</p>	<p>Maintenance Label + Barcode Label</p> <p>Position on the Product (Fig. 1 + 2)</p> <p>The maintenance label documents who did the regularly scheduled maintenance work and when.</p> <ul style="list-style-type: none"> - Mattress (Fig. 1) - Suspension framework (Fig. 2)

Remarks

General Remarks

- Read the Instructions for Use carefully before first use. If you have difficulty reading, please have them read to you. If the copy of the Instructions for Use is missing, you can download one from www.thomashilfen.de. The contents of the digital Instructions for Use can be enlarged.
- Check the packaging for any transport damage and to be sure the product is in flawless condition.
- Changes to the product are only permitted within the scope of the intended adjustment options.
- Cleaning and daily care can be provided by relatives and caregivers. Maintenance and inspection are done by the distributor.
- All labels have to remain on the product and should not be tampered with.
Exception: The maintenance label has to be updated and completed after each maintenance.

Safety Instructions

- Please do not overload the product and note the maximum user weight and the maximum load capacity (see technical data).
- The product should only be used by one user at a time. No other individuals should be placed in or on the product.
- Heavy perspiration, high humidity, and extreme changes in temperature can have an impact on the foam in the mattress. Pay attention to changes such as sagging in the mattress.
- If there is a serious and unexpected accident with the product, please inform our company immediately and the responsible authorities.

Intended Purpose

Intended Use / Application Areas

The ThevoSleeping Star Senso micro-stimulation system is a class 1 medical device.

The system can be used at home and in professional settings. It consists of a mattress and a system of suspensions with the MiS micro-stimulation effect that promotes natural movement and the user's body awareness. Special suspension blocks create a soft lying sensation. The special cut of the foam core also ensures a soft lying sensation and aids in air circulation to remove moisture. Overall, the system improves the lying and sleep quality.

The system replaces a standard mattress and goes into existing bed frames with the reclining surface and is also used in beds in care settings. The adjustment options of the existing bed frame are not limited by the micro-stimulation system. It is suitable for use in combination with digital assistance systems in care settings.

Any use that goes above and beyond this does not comply with the intended use.

Indication / Patient Target Group / Intended Users

The micro-stimulation system is suitable for active children and adolescents who sleep restlessly / wake up often and whose concentration and performance are limited during waking times as in the case of autism or ADHD.

This can definitely be used if there are illness-related difficulties falling asleep and sleep disturbances and impaired body awareness during sleep.

Risks Associated with Use and Contraindications

If used properly, risks can be eliminated during use. There are no known contraindications.

Product Description

Scope of Delivery

Your product is delivered with the standard equipment and comes preassembled. A copy of the Instructions for Use is enclosed with the product.

Standard Equipment

Mattress core; suspension framework; cover (jersey or incontinence)

Options

Upon request the product can be made in the following sizes:

- Width options from 60 cm to 100 cm (23.6" to 39.4")
- Length options from 120 cm to 200 cm (47.2" to 78.7")

The mattress is available for the following weight categories:

- ThevoSleeping Star Senso – up to max. 60 kg patient weight (132.3 lbs)

If the patient is in the upper end of a weight category (+/- 2 kg / 4.4 lbs), the next higher weight category should perhaps be chosen. We recommend a Thevo therapeutic mattress for adults.

If the patient's body has an atypical distribution of weight, the next weight category should perhaps be chosen.

Putting into Service

- Carefully remove the transport protection and the packaging.
- Keep the packaging material and small parts out of the reach of children to prevent swallowing or choking.
- The product can be used in any bed. Put the mattress with the red side up onto the slatted frame or your bed's supporting structure.
- The mattress can also be used without any slatted frame. In this case, there has to be a support surface of at least 13 cm / 5.1" in width on the long sides.
- Please note: For beds with side rails, side rail risers might be need.

Teaching the User

The user is taught how to use this device by trained staff who will reference these Instructions for Use.

Adjustment Options / Features

Accessories

Extra mattress cover (Fig. 3)

An extra mattress cover made of jersey or incontinence material is available if the original mattress cover has to be cleaned or discarded.

ThevoFloor – Floor bedding system (Fig. 4)

ThevoFloor is suitable for use at home and in professional care settings and allows patients to lie low to the floor. This reduces the need for measures which infringe on the freedom of the user and prevents falls out of bed. The system is comprised of:

- two wall parts (foam with incontinence cover)
- connecting belt
- cover strip for the hook and loop fastener
- A mattress and a system of suspensions with the MiS micro-stimulation effect are used that promote natural movement and the user's body awareness.

The system can be designed to suit the user. A standard Thevo mattress measuring 140 x 200 cm fits. Other sizes are available on request.

ThevoFloor is used in the corners of rooms. The Thevo mattress selected has to be appropriate for the indication and weight of the patient in question.

Put both wall parts on the walls in the corner of the room. Make sure that the hook and loop strip on the front of the short, top-wall part faces the hook and loop strips on the front of the longer side wall part (Fig. 5). Connect the wall parts firmly with the hook and loop strip. The cover strip can be put on the hook and loop fastener.

Pull the supplied connecting belt through the loops in the cover of the wall parts (Fig. 6).

Put the mattress red side up in the corner of the room. The top of the mattress should face the short part of the wall. Bring the connecting belt around the mattress. Adjust the length of the belt so that the wall parts are securely connected with the mattress (Fig. 7).

ThevoFloor wall parts	Top part	Side part
Length	140 cm / 55.1"	208 cm / 81.9"
Height	80 cm / 31.5"	80 cm / 31.5"
Width	8 cm / 3.1"	8 cm / 3.1"

Rescue belt system (Fig. 8)

With the rescue belt system bedridden patients can be removed quickly from danger. The rescue belt system is integrated into the cover of the mattress. The patient is secured by two hook and loop safety belts (at the chest and thigh). In an emergency, the hook and loop belts are pulled over the mattress, bedding, and patient and are connected. This way the patient and the bedding are secured quickly and easily to the mattress. The patient's arms are under the covers so that they are not injured during transport.

Using the pull straps at the top and foot of the bed the mattress, patient, and bedding are pulled out of the bed at a 90° angle. This way the patient can be evacuated horizontally and vertically from danger. The mattress slides very easily on all floor types. There are two lifting straps on each long side which allow emergency personnel to take the patient without transferring them from the bed.

The fabric of the belts is 100% polyester and is free of chemicals. This makes them non-irritant and also suitable for patients with allergies.

Maximum load capacity: up to 150 kg / 330.7 lbs

All users have to be familiar with how the product works. This way in an emergency the user can be removed quickly from the product.

Technical Data

ThevoSleeping Star Senso

Overall dimensions (w x l)	60 x 120 cm / 23.6 x 47.2"	70 x 140 cm / 27.6 x 55.1"	90 x 170 cm / 35.4 x 66.9"	90 x 200 cm / 35.4 x 78.7"	100 x 200 cm / 39.4 x 78.7"
Height	19 cm / 7.5"	19 cm / 7.5"	19 cm / 7.5"	19 cm / 7.5"	19 cm / 7.5"
Weight				18 kg / 39.7 lbs	20 kg / 44.1 lbs
max. user weight	60 kg / 132.3 lbs	60 kg / 132.3 lbs	60 kg / 132.3 lbs	60 kg / 132.3 lbs	60 kg / 132.3 lbs

Material Specifications

The open-pored and flame retardant cold foam provides moisture absorption and release and thermal regulation making the foam highly durable. The covers are flame retardant.

Cleaning, Maintenance, and Repair

Cleaning, Care Instructions, and Disinfection on a Daily Basis

Foam Mattress

The mattress can be autoclaved making it suitable for reuse.

Jersey Cover

The jersey cover is breathable and can be easily unzipped and removed for washing. Please also note the following care instructions.

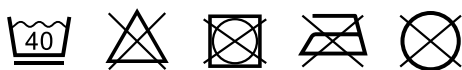
Incontinence Cover

The incontinence cover can be cleaned by wipe-down disinfection. The incontinence cover can be easily unzipped and removed for washing. Please also note the following care instructions.

Protective Cover for the Suspension framework

The black cover for the suspension framework can be easily unzipped and removed for washing. Please also note the following care instructions.

Suspension framework

Simple Cleaning: The bed base can be damp wiped with non-corrosive household cleaning agents or washed with a wet cloth. Make sure the bed base can dry well after. Cleaning with an alcohol-based disinfectant is also possible.

Disinfection

The mattress and suspension framework can be cleaned by steam vacuuming. This has to be done to disinfect and before each reuse.

We recommend two disinfecting methods recognised by the Robert Koch Institute:

- the fractionated vacuum method (vacuum-steam vacuum method) up to a maximum temperature of 110° C.
- the steam flow method up to a temperature of 110° C.

Important note: Thermal treatment above 110° C leads to premature ageing and deforming of the plastic parts which is why it shouldn't be done. When using chemical disinfectants always follow the information from the manufacturer to ensure they are suitable for plastics.

Maintenance and inspection by the distributor

Preparation

A copy of the Instructions for Use has to be enclosed with the product. If one is not available, request a copy from the distributor. Familiarise yourself with how the product works. If you are not familiar with the product, read the Instructions for Use before testing it. Clean the product before starting to test it. Please note the care instructions in the Instructions for Use and product-specific test instructions. If you have further questions, please contact the staff member you dealt with or Thomashilfen directly.

Service and repairs to the product should only be done by qualified staff. It is the user's responsibility to ensure that the prescribed inspections / maintenance are done regularly and in a timely fashion. If damage is apparent, the user has to actively inform and have the damage fixed immediately by qualified staff.

Maintenance Schedule – Annual Inspection

	Area	Function		Damage		Replace	
1	Cover	👁 / ↻		👁 / ↻			
2	Foam core	👁 / ↻		👁 / ↻			
3	Suspension framework	👁 / ↻		👁 / ↻			
4	Wing suspensions	👁 / ↻		👁 / ↻			
5	Suspension blocks	👁 / ↻		👁 / ↻			
6	Labels			👁 / ↻			
done by:				on:			
👁 = Inspection ↻ = Reuse ✓ = Test completed							

Spare Parts

Only original spare parts should be used for repairs. A list of the spare parts is available at www.thomashilfen.de. To receive the right spare parts or for possible technical recalls, we normally require the serial number / UDI-PI and the exact UDI-DI.

Removal / installation of spare parts should only be done by qualified staff! In the event of returns, please make sure the product is securely packaged for transportation.

Storage

Store the product in a safe place when it is not being used. Do not lay or put any heavy objects on the product while it is being stored. Protect the product from weather like rain, intense sunlight, and frost. If the product is wet, wipe the metal and plastic parts dry and let the padding dry before storing the product.

Disposal

If the product is no longer usable due to its condition or it has reached the end of its service life, it can be recycled at the local waste disposal company. Please refer to the local guidelines for recycling and disposal.

EU Declaration of Conformity

We hereby declare as the party solely responsible that this product and any versions and accessories comply with the basic regulations. The current EU Declaration of Conformity can be downloaded at www.thomashilfen.de or requested from us directly.

Warranty

Thomashilfen offers you the two-year warranty as required by law for this product. The warranty period starts when the item is handed over. The warranty includes all claims which impact function. Not included is damage which is the result of improper use (e.g. overloading) and natural wear and tear.

Period of Use

The product has a period of use of up to 10 years if used correctly and serviced regularly and professionally as per the information from the manufacturer. The product can be used beyond this point if it is in safe and technically sound condition. Time spent in storage at the distributor does not count towards the period of use.

The theoretical period of use does not equate to a guaranteed service life. It depends on the intensity of use and always undergoes an individual assessment by the distributor.

The period of use does not apply to wear parts like covers, the mattress core and so on.

Service Life

The product has a service life of 10 years or 2 reuses (3 users) if the requirements stipulated for the period of use are met. The service life consists of the period of use and the storage times between reuses.

Reuse

The product is suitable to be reused. Before being passed on to the next user the product always has to be carefully cleaned and disinfected. The condition has to be checked by trained staff for any wear and damage. Any issues have to be remedied to ensure the continued safe use of the product. A copy of the Instructions for Use is a component of the product and has to be passed on to the new user.

Distributor

Astris PME Pty Ltd

Young & Canberra Office:
Astris PME Pty Ltd
5 McVeigh St
Young NSW 2594

Sydney Office:
Astris PME Pty Ltd
30 Prime Dr
Seven Hills NSW 2147

Melbourne Office:
Astris PME Pty Ltd
9 Katherine Dr
Ravenhall VIC 3023

Phone: +61 1300 131 884
info@astris-pme.com.au
www.astris-pme.com.au

Lighthouse Medical

#3 - Golden drive
Coquitlam, British Columbia
V3K6R2
Customer Service: 886.609.6827
orders@lighthousemedicaltd.com
www.lighthousemedicaltd.com

O'Neill Healthcare Ltd.

P. O. Box 7196
Dublin 13
Phone: +353 1-8326509
Fax: +353 1-8326491
sales@onhealthcare.ie
www.onhealthcare.ie

Morton & Perry Ltd

Healthcare Equipment Solutions
Unit 4, 39 Porana Rd
Wairau Valley
North Shore
Auckland
Phone: +64 9 444 7285
dmorton@mortonperry.co.nz
www.mortonperry.co.nz

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